

## A BIOLOGICAL WARFARE DETECTION DEVICE (BIOWARD I)

Koen Gerritse<sup>1</sup>, Bram van der Gaag<sup>2</sup>, Edwin Stigter<sup>2</sup>, Stefan van Hövell tot Westerflier<sup>2</sup>  
Hans Harteveld<sup>1</sup> and Ruud Busker<sup>1</sup>

1. TNO-Prins Maurits Laboratory, Research Group Analysis of Toxic and Explosive Substances P. O. Box 45, 2280 AA Rijswijk, The Netherlands
2. TNO-Nutrition and Food Research Institute, Sensor Technology Department P.O. Box 360, 3700 AJ Zeist, The Netherlands

### ABSTRACT

Using a bulk acoustic wave immunosensor device, Staphylococcus enterotoxin B and two simulant Biological warfare agents (bacteriophage MS-2 and bacteria *Erwinia herbicola*) could be detected in a direct fashion without the use of labelled antibodies. The detection device is incorporated into a ruggedized suit-case together with an aerosol sampler and sample handling system and can be used in the field for military monitoring purposes.

### INTRODUCTION

Within the area of biological warfare there is a lack of modern biochemical sensors capable of real-time monitoring of pathogenic micro-organisms and toxins produced by them. General requirements for this type of sensors are: high sensitivity, short response times and high selectivity resulting in low false alarm rates. We here describe the development of an automated biosensor monitoring device for the detection of Staphylococcus Enterotoxin B (SEB), one of the protein toxins produced by the bacterium *Staphylococcus aureus*, MS-2 bacteriophage which is a simulant for a BW virus and *Erwinia herbicola* as a simulant for BW bacteria. These three agents are regarded as model agents covering a sufficiently broad spectrum.

The total system consists of an aerosol sampler, sample handling system, detection system and a PC and LCD screen for data acquisition. All subparts, which are thermostated if necessary, are mounted in a box, which can be transported by a single person. For detection purposes three 20 MHz piezoelectric quartz crystal sensor devices are employed. Anti-BW "simulant" antibodies are coupled to the sensor surface using several interfaces. The assay is based on a direct immunological detection by anti-BW "simulant" antibodies as catching molecules.

### RESULTS

By using the acoustic wave immunobiosensor system under laboratory circumstances, solutions of SEB or BW simulants could be detected in a direct fashion without the use of labelled compounds. The detection limit of the sensor for SEB was determined at 10 ng/ml. The detection limits for the simulant agents MS-2 and *Erwinia herbicola* were determined at 5 [x] 10<sup>7</sup> pfu and 10<sup>7</sup> cfu respectively. The sensor response is proportional to the concentration SEB or BW simulant in the samples.

After the sensor surface was cleaned by acid oxidation, an amino layer was deposited to the surface by plasma polymerisation. Antibodies directed to SEB or BW simulant were covalently bound to the amino functionalized surface [1]. Using the preconstructed interface the sensor response time for MS2 phages is approximately 10 min (Figure 1). SEB was injected as a control agent, which did not yield any response. The presented immunochemical detection principle is directly applicable to other BW agents.

## Report Documentation Page

*Form Approved*  
*OMB No. 0704-0188*

Public reporting burden for the collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing the collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden, to Washington Headquarters Services, Directorate for Information Operations and Reports, 1215 Jefferson Davis Highway, Suite 1204, Arlington VA 22202-4302. Respondents should be aware that notwithstanding any other provision of law, no person shall be subject to a penalty for failing to comply with a collection of information if it does not display a currently valid OMB control number.

1. REPORT DATE <b>00 JAN 2002</b>	2. REPORT TYPE <b>N/A</b>	3. DATES COVERED <b>-</b>			
4. TITLE AND SUBTITLE <b>A Biological Warfare Detection Device (Bioward I)</b>		5a. CONTRACT NUMBER			
		5b. GRANT NUMBER			
		5c. PROGRAM ELEMENT NUMBER			
6. AUTHOR(S)		5d. PROJECT NUMBER			
		5e. TASK NUMBER			
		5f. WORK UNIT NUMBER			
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) <b>TNO-Prins Maurits Laboratory, Research Group Analysis of Toxic and Explosive Substances P. O. Box 45, 2280 AA Rijswijk, The Netherlands; TNO-Nutrition and Food Research Institute, Sensor Technology Department P.O. Box 360, 3700 AJ Zeist, The Netherlands</b>		8. PERFORMING ORGANIZATION REPORT NUMBER			
		10. SPONSOR/MONITOR'S ACRONYM(S)			
9. SPONSORING/MONITORING AGENCY NAME(S) AND ADDRESS(ES)		11. SPONSOR/MONITOR'S REPORT NUMBER(S)			
		12. DISTRIBUTION/AVAILABILITY STATEMENT <b>Approved for public release, distribution unlimited</b>			
13. SUPPLEMENTARY NOTES <b>This article is from ADA409494 Proceedings of the 2001 ECBC Scientific Conference on Chemical and Biological Defense Research, 6-8 March , Marriott's Hunt Valley Inn, Hunt Valley, MD., The original document contains color images.</b>					
14. ABSTRACT					
15. SUBJECT TERMS					
16. SECURITY CLASSIFICATION OF:			17. LIMITATION OF ABSTRACT	18. NUMBER OF PAGES	19a. NAME OF RESPONSIBLE PERSON
a. REPORT <b>unclassified</b>	b. ABSTRACT <b>unclassified</b>	c. THIS PAGE <b>unclassified</b>	<b>UU</b>	<b>2</b>	

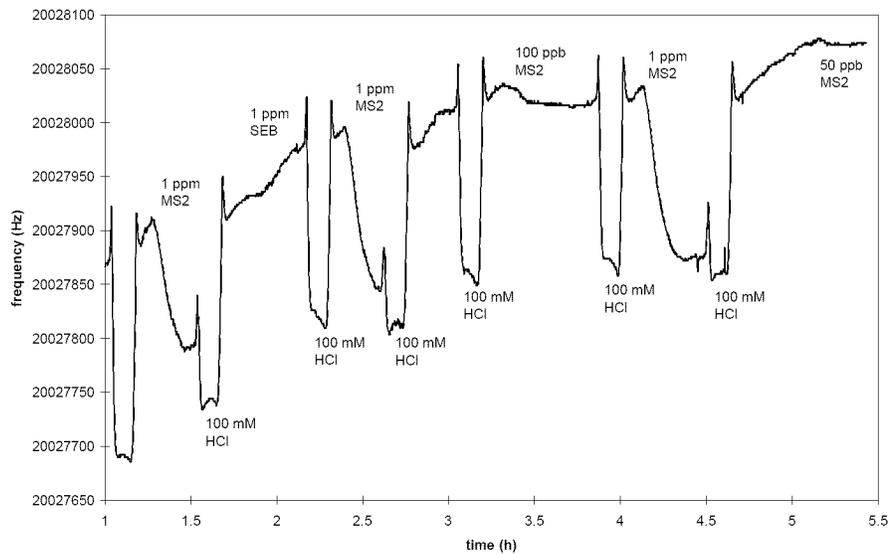


Figure 1. Sensor response of plasmapolymerised BAW crystals coated with anti-MS2 antibodies on injections of 1  $\mu\text{g}/\text{ml}$  MS2 in HBS and 1  $\mu\text{g}/\text{ml}$  SEB (control of specificity); 100  $\mu\text{l}$  sample injection in carrier flow rate of 25  $\mu\text{l}/\text{min}$ .

#### ACKNOWLEDGMENT

This work was performed in co-operation with GIAT Industries (France) and funded by the French and Dutch Ministries of Defence in the framework of Euclid RTP 13.7

#### REFERENCES

1. R. Nakamura, H. Muguruma, K. Ikebukuro, S. Sasaki, R. Nagata, I. Karube and H. Pedersen; A plasma polymerized film for surface plasmon resonance immunosensing; *Anal. Chem.* 69, 1997, 4649-4652.